

GOL SITE

GPP-100

Automatic Specific Protein Analyzer



FEATURES

GOLSITE



All in one
Cartridge



One Step




Reliable
Result



Compact

- Classical Immunoassay Technology: Nephelometry
- All-in-one cartridge for single test
- Auto ID of assay QR code on reagent cartridge
- Small sample volume and wide test menu
- Quantitative measurement and maintenance free
- Independent calibration

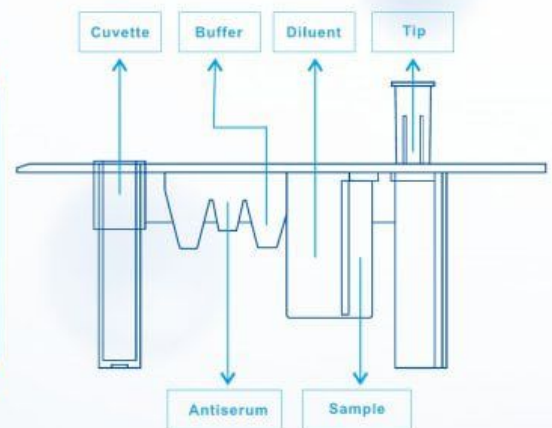

OPERATION

ONE STEP ONLY!



Add sample and insert
the cartridge, then just
wait for your test result!



TEST MENU




Assay	Linear Range	Sample Type	Category
HbA1c	2.8% - 15.5%	Blood	Diabetes
IgA	0.23 - 5.10g/L	Blood/Plasma/Serum	Immune Disease
IgG	1.8 - 26.0 g/L	Blood/Plasma/Serum	
IgM	0.18 - 4.60 g/L	Plasma/Serum	
C3	0.19 - 4.70 g/L	Plasma/Serum	
C4	0.06 - 2.20 g/L	Plasma/Serum	
KAP	0.88 - 12.20 g/L	Serum	
LAM	0.65 - 7.10 g/L	Serum	
ASO	41 - 607 IU/mL	Serum	
RF	20 - 306 IU/mL	Serum	Rheumatoid Disease
CCP	15 - 200 IU/mL	Serum	Cardiovascular Disease
D-Dimer	0.28 - 15.70mg/L	Plasma	
UsCRP	0.22 - 19.00mg/L	Plasma/Serum	
LP(a)	60 - 1200 mg/L	Serum	
FR-CRP	0.5 - 380.0 mg/L	Blood/plasma/Serum	Inflammation
CRP	1.6 - 152.0 mg/L	Plasma/Serum	
SAA	3 - 200 mg/L	Blood/plasma/Serum	
TRF	0.25 - 6.10 g/L	Serum	
Urine IgG	4 - 100 mg/L	Urine	
Urine TRF	1 - 80 mg/L	Urine	
A1M	5 - 200 mg/L	Urine	
mALB	9.5 - 230.0 mg/L	Urine	
CysC	0.4 - 7.6 mg/L	Plasma/Serum	
BMG	0.4 - 13.0 mg/L	Plasma/Serum/Urine	


More assays are to be expected

CLINICAL PERFORMANCE

GPP-100 VS Biorad D10

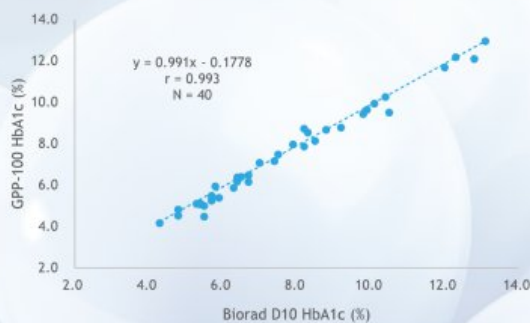

GPP-100 VS Beckman IMMAGE

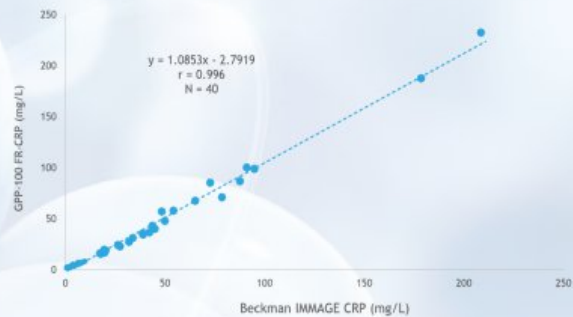

Excellent correlation result compared with well-established brand



SPECIFICATIONS

Measuring Time	3 to 8 minutes	Display	Color touch screen
Sample Volume	3 μ L to 300 μ L	Printer	Built-in thermal printer
Pack Size	20 T/Kit	Net Weight	7.1 kg
Memory Storage	Up to 100,000 test results	Size	404 mm*161 mm*397 mm

CONTACT US

Address:

No. 103C, 503C & 504D, Technology Building
 & No. 3A & 4A, Technology Building
 Annex, Zhaoshang Sub-District,
 Nanshan District, Shenzhen, China, 518067

Email: export@goldsite.com.cn

Tel: 0086 755 2689 0807

Website: en.goldsite.com.cn

Follow us on Twitter, YouTube, Facebook
 See our latest trend!



@ Goldsite Diagnostics